
AlphaNeb

IT > MANUALE ISTRUZIONI D'USO
 EN > INSTRUCTIONS FOR USE MANUAL
 DE > BEDIENUNGSANLEITUNG
 FR > MODE D'EMPLOI
 NL > GEBRUIKSAANWIJZING
 ES > MANUAL DE INSTRUCCIONES
 EL > ΕΓΧΕΙΡΙΔΙΟ ΟΔΗΓΙΩΝ ΧΡΗΣΗΣ
 ZH > 使用说明书
 < AR كتيب الاستخدام

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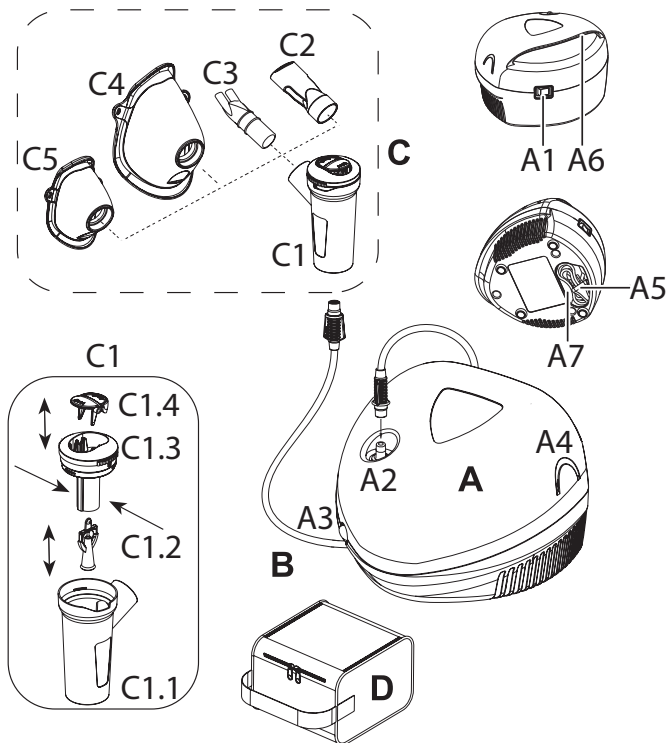Aerosol therapy device

ENGLISH

Congratulations on your purchase and thank you for choosing our unit. Our goal is to fully satisfy our customers by offering them cutting-edge products for treating ailments of the respiratory system. **Please read these instructions carefully and keep them handy for future reference. Use this device only as described in this manual. This is a home medical device to nebulise and administer medication prescribed or recommended by your doctor upon assessing the patient's general conditions.**

Please visit www.flaem.it to see our entire range of Flaem products.

THIS DEVICE COMES EQUIPPED WITH:

A - Aerosol Therapy Device

- A1 - Switch
- A2 - Air outlet port
- A3 - Air filter
- A4 - Nebulizer port
- A5 - Power cord
- A6 - Carrying handle
- A7 - Cable compartment porto

C - Accessories

- C1 - RF7 Dual Speed Nebulizer
 - C1.1 - Bottom piece
 - C1.2 - Nozzle
 - C1.3 - Top piece
 - C1.4 - Speed selector
- C2 - Mouthpiece
- C3 - Pediatric nose piece
- C4 - SoftTouch Adult mask
- C5 - SoftTouch Pediatric mask

B - Connection tube (main unit / nebulizer)

D - Convenient, roomy carrying bag

IMPORTANT SAFEGUARDS

- This device is also intended for direct use by the patient.
- Before using the product for the first time, and periodically during its lifetime, check the integrity of the device structure and of the power cable to make sure there is no damage. In the event of damage, do not plug in the cable and immediately take the product to an authorised FLAEM service centre or to your trusted dealer.
- Should your device fail to provide the expected performance, contact the authorised service centre for clarifications.
- The expected medical life of the accessories is 1 year. It is, however, advisable to replace the nebuliser cup every 6 months in the event of intense use (or earlier if the cup is obstructed) to always guarantee maximum therapeutic efficacy.
- Children and people who are not self-sufficient must use the device under the close supervision of an adult who has read this manual.
- Some parts of the device are small enough to be swallowed by children; therefore, keep the device out of the reach of children.
- Do not use the supplied tubing and cables for anything other than their intended use. These parts could cause a strangling hazard: pay close attention to children and persons with particular difficulties as they are often unable to accurately evaluate danger.
- The apparatus is unsuitable for use in presence of flammable anaesthetic mixture with air, oxygen or nitrous oxide.
- Always keep the power supply cable away from hot surfaces.
- Keep the power cable away from animals (for example, rodents) which could damage the insulation.
- Do not handle the device with wet hands. Do not use the device in damp environments (for example, while taking a bath or shower). Do not immerse the device in water; in the event of immersion immediately disconnect the plug. Do not remove or touch the

immersed device; unplug the power cable first. Immediately bring the device to an authorised FLAEM service centre or to your trusted dealer.

- Use the device only in dust-free conditions, otherwise treatment could be compromised.
- Do not wash the device under running water or by immersion and keep it safe from being sprayed by water or other liquids.
- Do not expose the device to particularly extreme temperatures.
- Do not place the device near sources of heat, in direct sunlight or in excessively hot rooms.
- Do not obstruct or put objects into the filter or its related housing in the device.
- Never obstruct the air vents located on both sides of the device.
- Always use it on a rigid surface that is clear of obstacles.
- Make sure there is no material obstructing the air vents before each use.
- Do not put any objects in the air vents.
- Repairs, including the replacement of the supply cord, are to be carried out by FLAEM authorised personnel only, by complying with the information provided by the manufacturer.
- The average expected duration for the compressor series are: F400: 400 hours, F700: 700 hours, F1000: 1000 hours, F2000: 2000 hours.
- **WARNING:** Do not modify this device without authorisation from the manufacturer.
- The Manufacturer, the Vendor and the Importer shall be held responsible for safety, reliability and performance only if: a) the device is used in compliance with the instructions for use b) the wiring where the device is being used is in compliance with safety regulations and current laws.
- Interactions: the materials used in contact with medication have been tested with a vast range of medications. However, in view of the variety and continuous evolution of pharmaceuticals, interactions cannot be ruled out. We recommend using the medication as soon as possible once it has been opened and preventing prolonged exposure in the nebuliser cup.
- The manufacturer must be contacted about any problems and/or unexpected events concerning operation and for any clarifications on use, maintenance/cleaning.
- Interactions: The materials used in the medical device are biocompatible in accordance with the provisions of Directive 93/42 EC and subsequent amendments. However, the possibility of occurrence of allergic reactions cannot be entirely excluded.
- The amount of time required between preservation conditions and use is approximately 2 hours.

OPERATING INSTRUCTIONS

Before each use, clean hands thoroughly and clean the device as described in the section on “CLEANING, SANITIZATION AND DISINFECTION”. During use, it is advisable to protect yourself from any dripping. It is recommended that each person use their own nebulizer cup and accessories to prevent risk of infection due to contamination. This device is suitable for the administration of medical substances and not, for which the administration via aerosol is foreseen; these substances are to be in any case prescribed by the Doctor. In case of too thick substances, the dilution with a suitable physiological solution could be needed, according to the medical prescription.

1. Plug the power cord (A5) into a power socket corresponding to the voltage of the unit. This must be positioned so that it is not difficult to make the disconnection from the mains.
2. Insert the nozzle (C1.2) in the upper part (C1.3) pressing as shown by the 2 arrows in the “Connection diagram” in point C1. Insert the Speed selector (C1.4) in the upper part (C1.3) as shown in the “Connection diagram” in point C1. Pour the medication prescribed by the doctor into the lower part (C1.1). Close the nebuliser by turning the upper part (C1.3)

clockwise.


3. Connect accessories as indicated in the "Connection diagram" on the cover.
4. Sit comfortably holding the nebulizer in your hand, place the mouthpiece onto your mouth or alternatively use the nose piece or mask. Should you opt for the mask accessory, place it on your face as shown in the picture (with or without using the elasticated strap).
5. Turn on the device using the switch (A1) and then breathe in and exhale deeply.
6. When the treatment is finished, shut off the device and unplug the power cord.


WARNING: If after the therapy session moisture accumulates in the tube (B), remove the tube from the nebulizer and dry it using the device's compressor fan; this operation will prevent mold from forming inside the tube.

HOW TO USE THE "RF7 DUAL SPEED" NEBULIZER WITH THE SPEED SELECTOR

Professional and fast, this device is suitable for administering all types of medications, including more expensive ones, even in patients with chronic diseases. The geometry of the internal lines of the RF7 Dual Speed nebulizer ensures the ideal granulometry for active treatment all the way down to the lower respiratory tract.

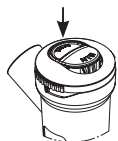
To speed up the inhalation therapy, move the speed selector button (C1.4) by pressing on MAX with your finger.

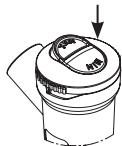
To increase the effectiveness of the inhalation therapy, press with your finger on the wording MIN of the speed selector (C1.4). In this position, the speed selector acts as a valve and allows nebulizing the optimal amount of drug to the lower respiratory tract, reducing its dispersion in the environment.


SOFTTOUCH FACE MASKS

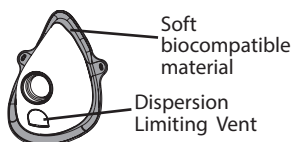SoftTouch masks have an outer edge made of **soft biocompatible material** that ensures excellent adherence to the face, and is also equipped with an **innovative Dispersion Limiting Device**. These distinctive elements that distinguish it allow greater sedimentation of medication in the patient and also **limit dispersion**.

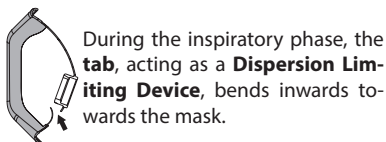During the inspiratory phase, the **tab**, acting as a **Dispersion Limiting Device**, bends inwards towards the mask.

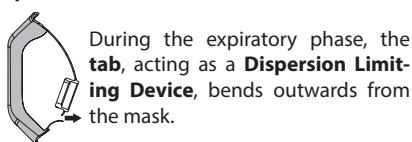During the expiratory phase, the **tab**, acting as a **Dispersion Limiting Device**, bends outwards from the mask.

CLEANING SANITISATION DISINFECTION

Switch off the device before any cleaning procedure and unplug the power cable from the socket.

DEVICE AND TUBING EXTERIOR

Use only a damp cloth with antibacterial soap (non-abrasive and with no solvents of any sort).

ACCESSORIES

Open the nebuliser by turning the upper part (C1.3) anticlockwise, remove the nozzle (C1.2) from the upper part (C1.3) by pressing as shown by the 2 arrows in the "Connection diagram" in point C1.

Then proceed according to the following instructions.

SANITISATION

Before and after each use, sanitise the nebuliser cup and the accessories, choosing one of the methods described below.

method A: Sanitise accessories **C1.1-C1.2-C1.3-C1.4-C2-C3-C4-C5** under potable hot water (approximately 40°C) with a gentle, non abrasive dish detergent.

method B: Sanitise accessories **C1.1-C1.2-C1.3-C1.4-C2-C3-C4-C5** in the dishwasher with a hot cycle.

method C: Sanitise accessories **C1.1-C1.2-C1.3-C1.4-C2-C3-C4-C5** by immersing them in a solution of 50% water and 50% white vinegar, then rinse thoroughly under potable hot water (approximately 40 °C).

If you want to also perform the cleaning for DISINFECTION, jump to the DISINFECTION paragraph.

After having sanitised the accessories, shake them vigorously and place them on a paper towel. Alternatively, dry them with a jet of hot air (for example, a hair dryer).

DISINFECTION

After sanitising the nebuliser cup and the accessories, disinfect them choosing one of the methods described below.

method A: Accessories **C1.1-C1.2-C1.3-C1.4-C2-C3-C4-C5** can be disinfected.

The disinfectant must be an electrolytic chloroxidizer (active principle: sodium hypochlorite) specific for disinfecting, which is available in any pharmacy.

Implementation:

- Fill a container big enough to hold all of the parts to disinfect with a solution of potable water and disinfectant, according to the proportions indicated on the packaging of the disinfectant.

- Completely immerse each part in the solution, taking care to avoid the formation of air bubbles on the parts. Leave the parts immersed for the amount of time indicated on the packaging of the disinfectant associated with the concentration chosen for the solution.
- Remove the disinfected parts and rinse abundantly with warm potable water.
- Dispose of the solution following the instructions provided by the disinfectant manufacturer.

method B: Sanitise the accessories **C1.1-C1.2-C1.3-C1.4-C2-C3-C4-C5** by boiling them in water for 10 minutes; use demineralised or distilled water to prevent calcium deposits.

method C: Sanitise the accessories **C1.1-C1.2-C1.3-C1.4-C2-C3-C4-C5** with a hot steam steriliser for baby-bottle (not the microwave type). Perform the process faithfully following the instructions of the steriliser. To ensure that the disinfection is effective, choose a steriliser with an operating cycle of at least 6 minutes.

After having disinfected the accessories, shake them vigorously and place them on a paper towel. Alternatively, dry them with a jet of hot air (for example, from a hair dryer).


At the end of each use store the device complete with accessories in a dry place away from dust.

REPLACING THE FILTER

The device has a suction filter (A3) that should be replaced when it is dirty or changes color. Do not wash or reuse the same filter. Regular replacement of the filter is required to help and ensure correct performance of the compressor. Check the filter on a regular basis. Do not replace the filter during use. Contact your local dealer or authorized service center for spare filters.

To replace the filter:

Insert a flathead screwdriver between the edge of the filter and the body. Lift the filter and remove it by turning and pulling it. The filter has been designed so that it is always fixed in its seat.

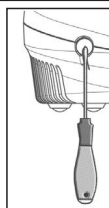AVAILABLE SPARE PARTS

Description

- RF7 Dual Speed Nebulizer Set composed of: nebulizer and mouthpiece
- Mask Set composed of: adult mask, pediatric mask and elastic band
- 1 m connection tube
- Spare Air filter kit per unit (2 pcs)
- Nose piece Set composed of: adult nose piece and pediatric nose piece

Code

ACO461P
ACO462P
ACO35
ACO439P
ACO463P

TROUBLESHOOTING

Switch off the device before any procedure and unplug the power cable from the socket.

PROBLEM	CAUSE	SOLUTION
The device does not work	The power cable has not been correctly inserted into the power socket	Correctly insert the power cable in sockets.
The device does not nebulize or nebulizes insufficiently.	The medication has not been inserted in the nebulizer	Pour the right amount of medication in the nebulizer
	The nebulizer has not been properly fitted.	Disassemble the nebulizer and reassemble it correctly as shown in the connection diagram on the cover.
	The nebulizer nozzle is clogged	Disassemble the nebulizer, remove the nozzle and perform cleaning operations Failure to remove medication deposits from the nebulizer affects its efficiency and operation. Strictly comply with the instructions contained in the “CLEANING, SANITISATION, DISINFECTION” chapter
	The air pipe is not correctly connected to the device	Ensure the device’s accessories are properly connected to the air intake (see connection diagram on the cover).
	The air pipe is bent, damaged or twisted	Unwind the pipe and ensure it is not crushed or punctured. If necessary replace it.
	The air filter is dirty	Replace the filter
The apparatus is noisier than usual	The Filter is not correctly inserted in its housing	Fully and properly insert the Filter in its housing

If after verifying the above mentioned conditions the device should not operate properly, contact your trusted retailer or the nearest authorized FLAEM service center.

SYMBOLS

Certification TÜV

CE Medical Marking ref. Dir. 93/42
EEC and subsequent amendments

Class II device

Important: check the operating
instructions

"OFF" for part of equipment

"ON" for part of equipment

Complies with: European standard
EN 10993-1 "Biological Evaluation
of Medical Devices" and European
Directive 93/42/EEC "Medical
Devices". Phthalates free. In
conformity with Reg. (EC) no.
1907/2006

Minimum and maximum room
temperature

Minimum and maximum
atmospheric pressure

Serial number of device

Manufacturer

Type BF applied part

Alternating current

Enclosure protection rating: IP21.
(Protected against solid bodies over
12 mm. Protected against access with
a finger;
Protected against vertically falling water
drops).

Minimum and maximum air moisture

ELECTROMAGNETIC COMPATIBILITY

This device was designed to satisfy the currently required requisites for electromagnetic compatibility (EN 60601-1-2). Electro-medical devices require particular care during installation and use relative to EMC requirements. Users are therefore requested to install and/or use these devices following the manufacturer's specifications. Risk of potential electromagnetic interference with other devices. Mobile or portable RF radio and telecommunication devices (mobile phones or wireless connections) may interfere with the operation of medical devices. For more information, visit the website www.flaemnuova.it. The appliance may be susceptible to electromagnetic interference in the presence of other devices used for specific diagnoses or treatments. Flaem reserves the right to make technical and functional changes to the product without prior notice.

TECHNICAL FEATURES

Mod. P0707EM F400

Voltage:	230V~ 50Hz 130VA	115V ~ 60Hz	220V ~ 60Hz
Safety certifications:			

Max pressure:	1.8 ± 0.3 bar
Compressor air output:	9 l/min approx
Sound level (at 1 m):	54 dB (A) approx
Operation:	Continuous use

Operating conditions:

Temperature:	min 10°C; max 40°C
RH Air humidity:	min 10%; max 95%
Atmospheric pressure:	min 69KPa; max 106KPa

Storage conditions:

Temperature:	min -25°C; max 70°C
RH Air humidity:	min 10%; max 95%
Atmospheric pressure:	min 69KPa; max 106KPa

Dimensions (W)x(D)x(H):	18x19x13 cm
Weight:	1.300 Kg

APPLIED PARTS

Type BF applied parts are: patient accessories (C2, C3, C4, C5)

RF7 Dual Speed nebulizer

Medication minimum capacity:	2 ml
Medication maximum capacity:	8 ml
Operating pressure (with neb.):	0.65 bar

	Speed selector C1.4	
	Max	Min
⁽¹⁾ Delivery:	0.55 ml/min approx.	0.25 ml/min approx.
⁽²⁾ MMAD:	4.58 µm	3.78 µm
⁽²⁾ Breathable fractions < 5 µm (FPF):	54.4%	63%

⁽¹⁾ data detected according to Flaem I29-P07.5 internal procedure

⁽²⁾ In vitro testing performed by TÜV Rheinland Italia S.r.l. in collaboration with the University of Parma and in compliance with the EN 13544-1:2007 + A1 European Standard for aerosol therapy devices Standard. More details are available on request.

DISPOSAL OF DEVICE:

In conformity with Directive 2012/19/EC, the symbol shown on the device to be disposed of indicates that it is considered as waste and is therefore subject to "sorted waste collection". The user must therefore take (or have taken) the above waste to a pre-sorted waste collection centre set up by the local authorities, or else give it back to the dealer when purchasing a new appliance of the same type. Pre-sorted waste collection and the subsequent treatment, recovery and disposal operations favor the production of appliances made of recycled materials and limit the negative effects of any incorrect waste management on the environment and public health. The unlawful disposal of the product by the user could result in administrative fines as provided by the laws transposing Directive 2012/19/EC of the European member state or of the country in which the product is disposed of.

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© 2020 FLAEM NUOVA®
FLAEM NUOVA S.p.A. non risponde di ulteriori estensioni del periodo di garanzia assicurate da terzi.
FLAEM NUOVA S.p.A.